


# W+File

One File Endo

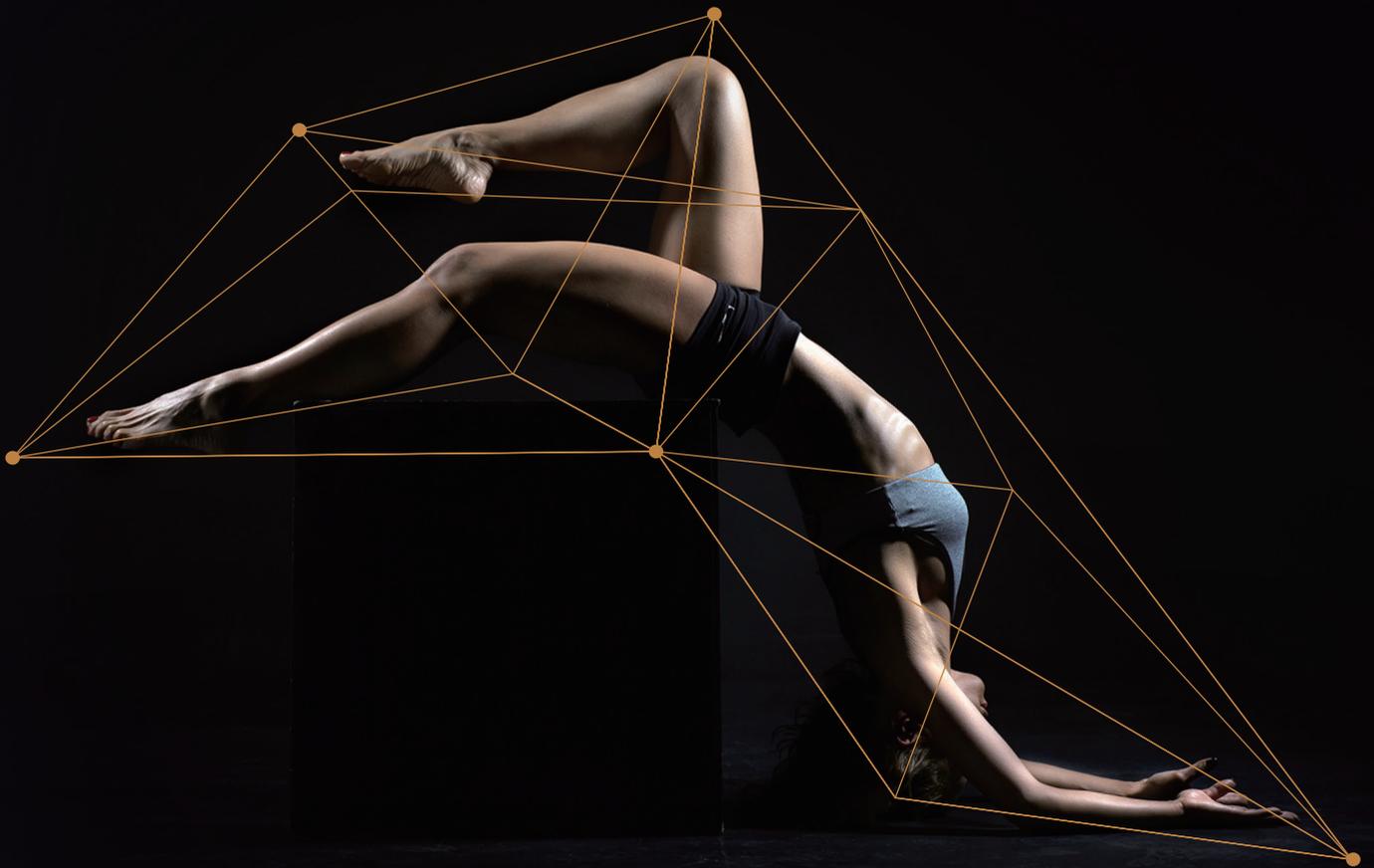

**NIC**<sup>®</sup>

For Better Endo

# W+File

One File Endo




Radiopaque stopper



Precise depth marking




Reciprocating

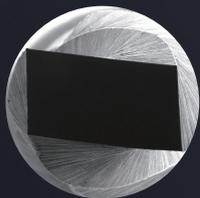

Parallelogram



Safe guiding tip

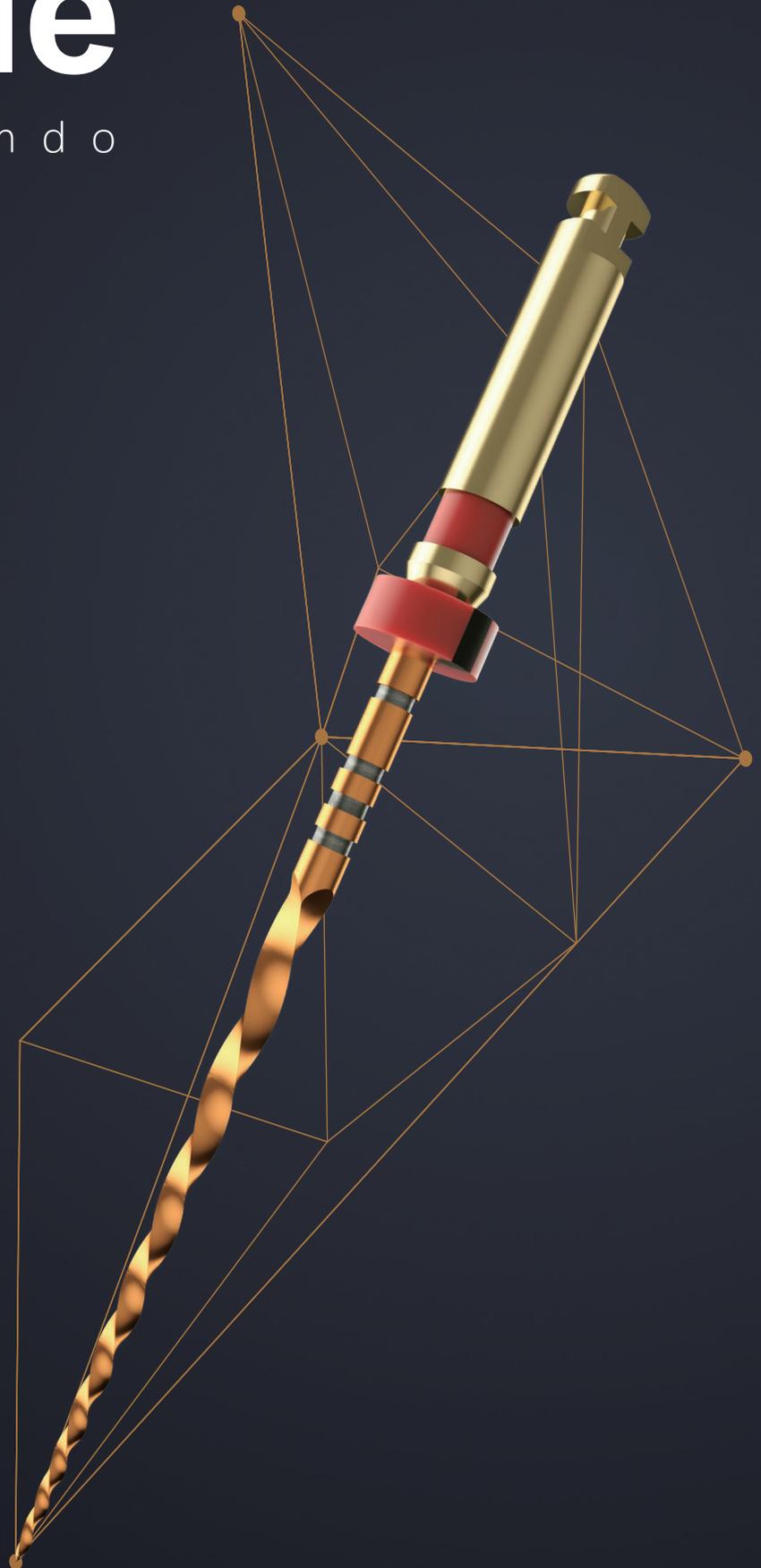

# Features

- MaxTech 3.0 NiTi alloy
- Parallelogram cross section
- Variable regressive taper design to have accurate shaping
- Strong cutting ability shorten cutting time
- More space for debris removal
- Progressive pitch reduce screw-in
- One file system



W+File	Small	Primary	Medium	Large
Size	20/.07	25/.07	35/.06	45/.05
Speed(RPM)	Reciprocating	Reciprocating	Reciprocating	Reciprocating
Torque(N.CM)	Reciprocating	Reciprocating	Reciprocating	Reciprocating
GP	Small	Primary	Medium	Large

★ If Speed is adjustable under reciprocating motion, V file is recommended to set at:300-350RPM.

# Sequence

- Scout the canal with hand use file
- Create glide path by any path function file
- Choose proper size shaping files
- Shaping the canal with the chose W+File until it reaches working length
- After shaping and finishing is done, then choose the matched GP points for obturation.

SSK #10



T+File



Primary



◆ Small 20/.07 for narrow canals

◆ Primary 25/.07 for most cases

◆ Medium 35/.06 & Large 45/.05  
for wide canals



Small  
20/.07



Medium  
35/.06



Large  
45/.05

# References

Size	21 mm	25 mm	31 mm
20	AFW 021 020	AFW 025 020	AFW 031 020
25	AFW 021 025	AFW 025 025	AFW 031 025
35	AFW 021 035	AFW 025 035	AFW 031 035
45	AFW 021 045	AFW 025 045	AFW 031 045
20/25/35/45	AFW 021 A01	AFW 025 A01	AFW 031 A01

# Clinic Cases



# Endo Motor & Obturation



Model	Ref. No
Small	PPW 029 020
Primary	PPW 029 025
Medium	PPW 029 035
Large	PPW 029 045
Assorted	PPW 029 B01



Model	Ref. No
Small	GPW 029 020
Primary	GPW 029 025
Medium	GPW 029 035
Large	GPW 029 045
Assorted	GPW 029 B01

Premium manufacturer  
since 2002

SUPERLINE



SHENZHEN SUPERLINE TECHNOLOGY CO., LTD.

ISO 13485:2016 CE 197 FDA

[www.supline.com](http://www.supline.com)